

FULL PROTECTION

APRON FEATURES

+ Adjustable Buckle Set
Adjustable buckle closures ensure a secure fit.

+ Non-slip Backing (skirt interior)
Reduces slippage and improves fit

+ Extended Front Panels
Front panels are extended to 3/4 overlap coverage where LE protection is doubled

Shoulder Padding
Premium foam padding in the shoulders provides extra comfort.

Apron Pocket

Velcro and Buckle Closures

Hanging Tabs on Skirt

ADDITIONAL CUSTOMIZATION OPTIONS


Thyroid Collar Wings Extra Pocket Suspenders (skirt)
Embroidery Logo Sleeve

MODEL: [AP-VQ]_(vest only) [AP-SKIRT]_(skirt only) [AP-VSQ]_(vest and skirt combo)

Be prepared for prolonged procedures with the Vest Skirt apron. The two piece design shifts 70% of the weight directly to the hips, reducing back fatigue. The Essential apron configuration **features extended front panels to achieve a 3/4 overlap (where protection is doubled) which allows for a lighter weight apron without compromising protection.** The skirt comes standard with a premium non-slip fabric in the back. **Essential aprons are only available in standard sizes.** For a made-to-measure Vest Skirt, see the Balance Vest Skirt Apron.



Example of a 0.35mm LE Essential Vest Skirt Apron Lead Configuration

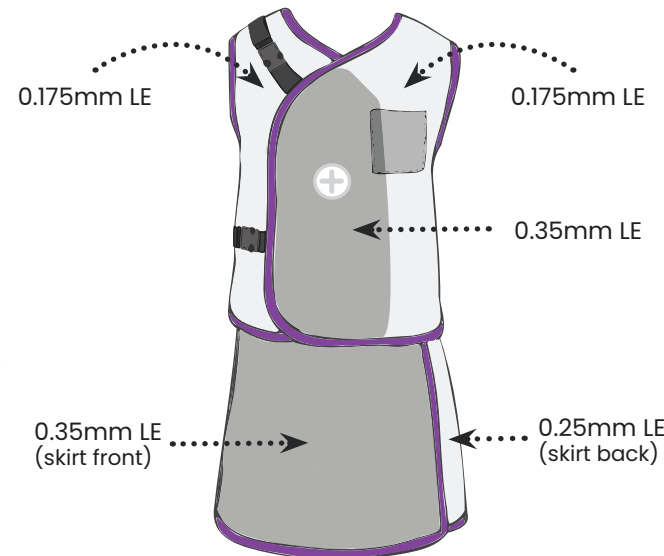

Core Material	Weight (0.25mm)	Thickness Options
Lead	3500 gsm	0.25, 0.50mm
Prolite	2850 gsm	0.25, 0.35, 0.50mm
Prolite IEC	3000 gsm	0.25, 0.35, 0.50mm
Premier	2600 gsm	0.25, 0.35, 0.50mm

For sizing and color options, please see our sizing and fabric color documents or view them online at protechmed.com.

*CE applies to Lead and Prolite IEC in 0.25, 0.35, and 0.50mm. Prolite IEC is limited to 130kVp diagnostic range.